
NAUDOJIMO INSTRUKCIJA

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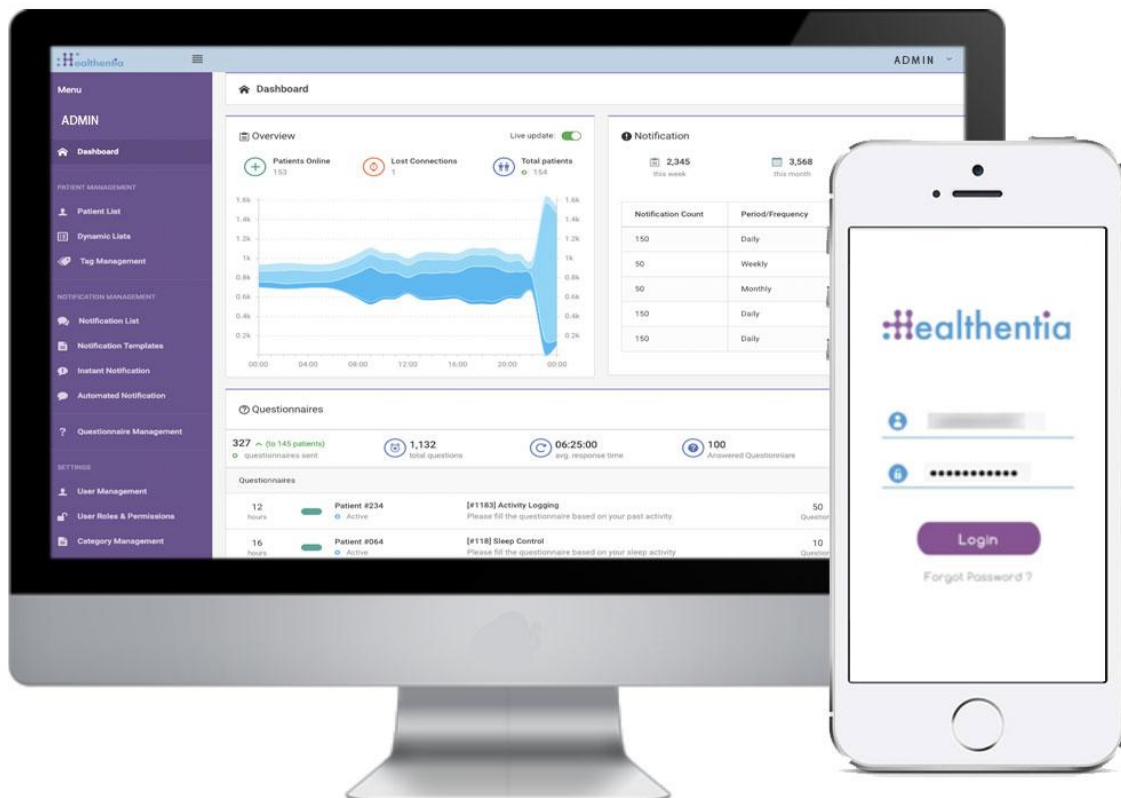Icon	Name	Organization	Codename	Disease	Patients	Status	Alerts	Actions
	Test Study-I	ORG_1	testirida		0	ONGOING	0	
	Interface	ORG_6	Interface	HIV Remote Patient ...	66	ONGOING	10	
	Symptom Tracker	ORG_1	COVID19		356	ONGOING	285	

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Account Settings

Name: * Therapeutic Area: Codename: *

Protocol: Duration In Months: Start Date: 2021/05/14 End Date:

Primary Language: * Add More Languages: Study Logo: 1

--Select-- Choose File No file chosen

Additional Services

☐ Video-Chat ☐ Multicenter ☐ Interactive Map

MOBILE CONFIGURATION

Study Information Page: Study Logo on App: 1

Leave empty for no Information Page Choose File No file chosen

Pop-up for Study's terms & Privacy

Short Term's Description: Terms of Use URL: 1

Leave empty for no Short Terms Description Leave empty for no Terms of Use

Privacy Policy URL: 1

Leave empty for no Privacy Policy

Available Consents

☒ General Terms ☒ Allowing Reporting ☐ Contacted by doctor in case of need

Additional App Widgets

☐ Treatment Reminder ☐ Two-Factor Authentication ☐ Liquid consumption ☐ Sleep Widget

Activity Widget (Available integrations)

☒ Fitbit Sync ☒ Garmin Sync ☒ Apple Health Sync ☒ Android Sync

SUBJECT DETAIL PAGE CONFIGURATION

☐ First/Last Name ☐ Fiscal Number ☐ Study Dates ☐ Severity

☐ Managed Fitbit ☐ Tracker Description ☐ Mobile Description ☐ Therapeutic Area

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3.1 INTERACTIVE MAP

At the configuration of the study, the Study Admin can select to enable the interactive map if there are questions regarding location. Portal users can select parameters through several filters and have an overview of subject outcomes at a geographical level, as presented in Figure 6.



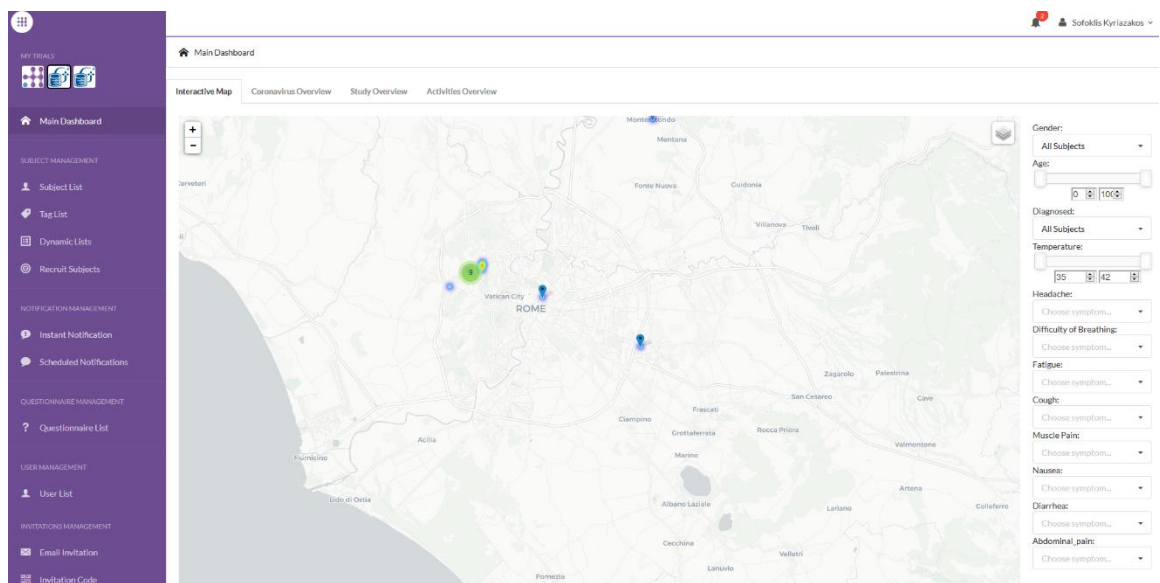
Figure 6: Interactive Map

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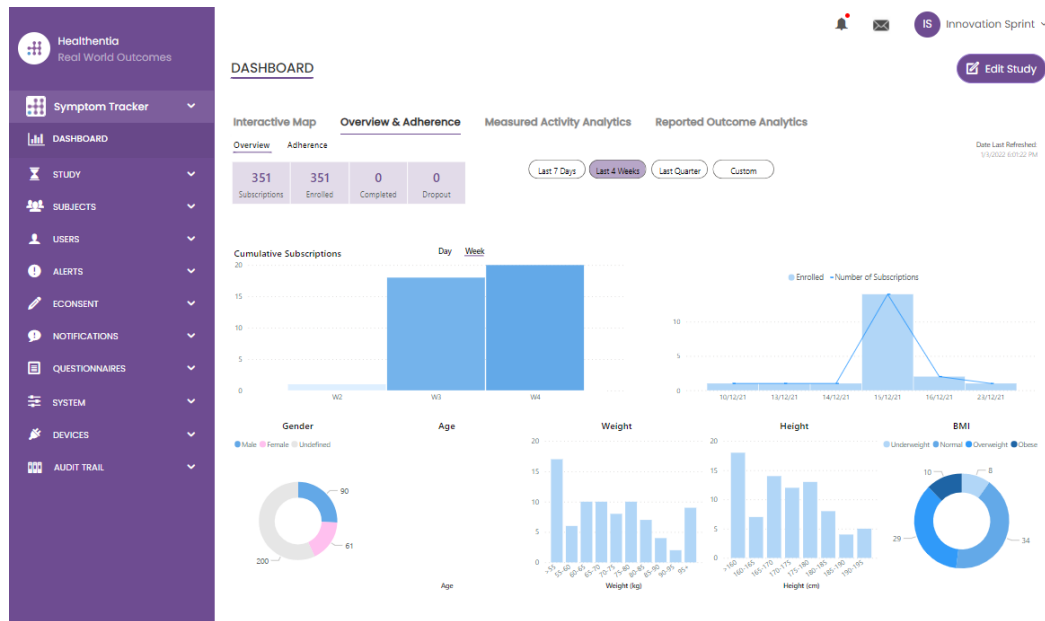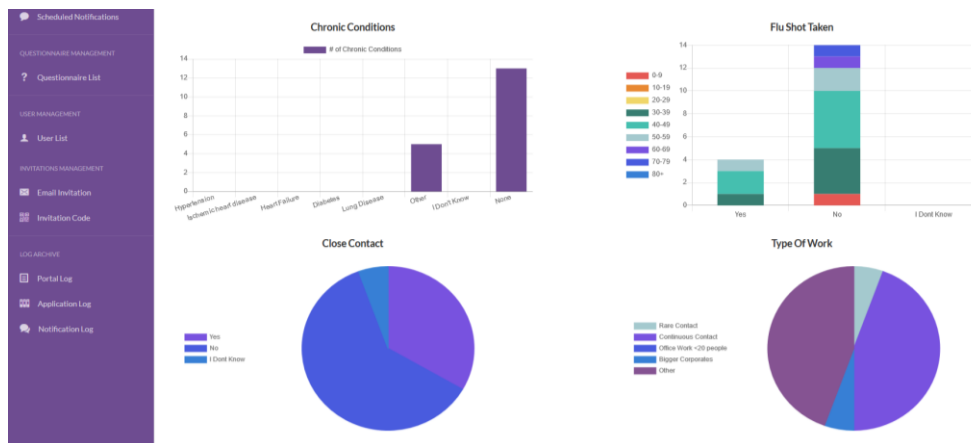The screenshot shows a 'Create Invitation Code' modal dialog box. It contains the following elements:

- Title:** Create Invitation Code
- Text:** Type a Codename or generate a random name by using the icon
- Form:**
 - A text input field labeled 'Code' with a random name icon to its right.
 - A numeric input field labeled 'Number of times the Code can be used:' with a 'Total times' placeholder.
 - A checkbox labeled 'Used together with an invited email'.
- Button:** A 'Create Code' button with a plus icon.

The background shows a sidebar with 'Healthlanta Real World Outcomes', 'Symptom Tracker', 'DASHBOARD', 'STUDY', 'SUBJECTS', 'Subject List', 'Tag List', 'Dynamic Lists', 'Custom Fields', 'Subject Invitations', 'Codes', and 'USERS'. The main area displays a table of 'Subject Invitations' with columns for 'Active', 'Created', 'Code', 'Number of times the Code can be used', and 'Status'.

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The screenshot shows the 'TAG MANAGEMENT' section of the Healthentia web portal. On the left is a purple sidebar with navigation options: Healthentia Real World Outcomes, Symptom Tracker, DASHBOARD, STUDY, SUBJECTS, Subject List, Tag List (selected), Dynamic Lists, Custom Fields, Subject Invitations, and Codes. The main area has a header with 'TAG MANAGEMENT' and a search bar. Below is a table with columns: Title, Created On, Type, Status, and Actions. The table contains six rows of tags: Low, Medium, High, Dev, and Not Synced, all with a status of 'ACTIVE'. At the bottom, there are pagination controls showing 5, 10, and All items, with the first page selected.

Title	Created On	Type	Status	Actions
Low	23/03/2020	Auto	ACTIVE	[Edit] [Delete]
Medium	23/03/2020	Auto	ACTIVE	[Edit] [Delete]
High	23/03/2020	Auto	ACTIVE	[Edit] [Delete]
Dev	01/04/2020	Custom	ACTIVE	[Edit] [Delete]
Not Synced	29/05/2020	Auto	ACTIVE	[Edit] [Delete]

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The screenshot shows the 'CREATE DYNAMIC ALERT' form in the Healthentia web portal. The sidebar is the same as in the previous screenshot. The main area has a header with 'CREATE DYNAMIC ALERT'. The form includes fields for: Name (with an asterisk), Status (dropdown menu set to 'Active'), Portal Alert Message (with a hint about using {shortid}), Push Notification Message, Email List (with a hint about separating emails), Email Message (with a hint about using {shortid}), API Endpoint to Call, API Alert Codename, and Tags to Apply. At the bottom right are 'Cancel' and 'Create' buttons.

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ALERTS

Filter by Tags: COVID19 POSITIVE, DOSE 1 COVIDVAX, DOSE 2 COVIDVAX, FEEDBACK, FITBIT TOKEN EXPIRED, HIGH COVID RISK, HIVSRQ SYMPTOM, LOW ADHERENCE, LOW COVID RISK, NOT SYNCED, PAZIENTE FRAGILE, PSYCHIATRIC REDDIT, PSYCHOLOGICAL SUPPORT

Show: 50

Triggered On (Browser Time)	Subject Id	Alert Message	Tag	Manager
21/04/2021, 14:21:27	20575866	Patient 20575866 needs psychological support	PSYCHOLOGICAL SUPPORT	Giulia Micheli
20/04/2021, 16:30:02	14455901	Subject 14455901 has «Questionario sulla Salute - EQ-5D-3L» as pending		Giulia Micheli
18/04/2021, 16:30:01	04588190	Subject 04588190 has «Questionario sulla Salute - EQ-5D-3L» as pending		Giulia Micheli
18/04/2021, 13:15:19	14455901	Patient 14455901 is identified as FRAIL	PAZIENTE FRAGILE	Giulia Micheli
18/04/2021, 13:12:09	14455901	Subject 14455901 has an HIVSRQ symptom alert	HIVSRQ SYMPTOM	Giulia Micheli
17/04/2021, 19:19:56	03806999	Patient 03806999 needs psychological support	PSYCHOLOGICAL SUPPORT	Giulia Micheli
17/04/2021, 16:30:02	03806999	Subject 03806999 has «Questionario sulla Salute - EQ-5D-3L» as pending		Subject answered the questionnaire.
15/04/2021, 16:30:01	15929649	Subject 15929649 has «Questionario sulla Salute - EQ-5D-3L» as pending		Giulia Micheli
10/04/2021, 09:10:41	11294295	Patient 11294295 needs psychological support	PSYCHOLOGICAL SUPPORT	Cristina Seguiti

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DYNAMIC LISTS

Drag a column header here to group by that column

Search...

Title	Description	Created On	Status	Actions
test	test	05/03/2021	ACTIVE	

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Healthentia
Real World Outcomes

Symptom Tracker

DASHBOARD
STUDY
SUBJECTS
USERS
ALERTS
ECONSENT
NOTIFICATIONS
QUESTIONNAIRES
SYSTEM
DEVICES
AUDIT TRAIL

DETAILS DYNAMIC LISTS

Title:
test
Status:
Active

Description:

Patient Criteria

Status:
Enrolled
Gender:
---Select---
Age:

Baseline Date:
Termination Date:
Tags:
Not Synced

Show:
10

Sr. No	Subject ID	Registration Date
No matching records found		

Showing 0 to 0 of 0 entries (filtered from 95 total entries)

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Description	Type	Content	Recipients	Sent On (User Time)	Sent On (Your Time)
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	02/02/2022, 16:24:58
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	02/02/2022, 14:21:49
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	31/01/2022, 16:28:38
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	31/01/2022, 16:27:03
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	28/01/2022, 20:46:50
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	28/01/2022, 12:33:32
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	27/01/2022, 15:43:10
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	27/01/2022, 14:39:19
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	21/01/2022, 15:03:14
<!--Instant Notification-->	Instant	View Content	View Recipients (1)	<!--Based on each user's time-->	18/01/2022, 11:54:26

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SCHEDULED NOTIFICATION

Drag a column header here to group by that column

Description	Type	Content	Status	Actions
Weekly COVID19 Self-Assessment	Weekly	View Content	INACTIVE	Edit Delete
Test	Monthly	View Content	INACTIVE	Edit Delete

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EDIT SCHEDULED NOTIFICATION

Description: Status:

Schedule Type: ☒ Regular ☐ By Baseline Date

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Subject	Title	Last Message	Sent On	Category	Reply
SBCYT	Information	Subject: Like arm pain	04/05/2021, 18:20:28	HEALTH	

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CONVERSATION WITH KONSTANTINA KOSTOPOULOU: ΔΩΝΒΓ

START TELECONSULTATION

Apr 22, 2021

hello 12:54 PM

Teleconsultation
Teleconsultation request was sent. 12:55 PM
Teleconsultation has been completed. 12:57 PM

Teleconsultation
Teleconsultation request was sent. 12:58 PM
Teleconsultation has been completed. 12:58 PM

Teleconsultation
Teleconsultation request was sent. 6:57 PM
Teleconsultation has been completed. 6:58 PM

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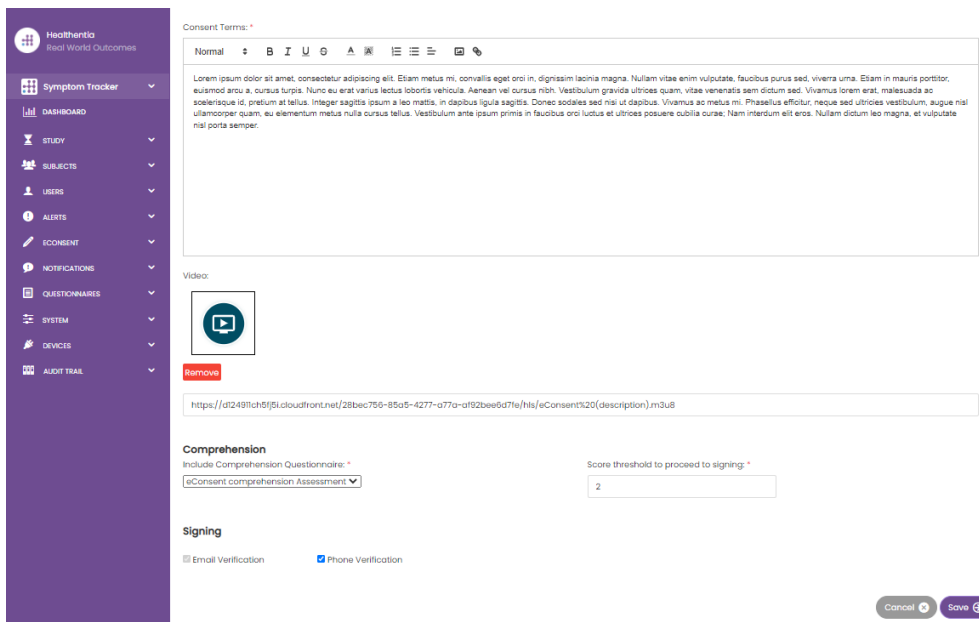
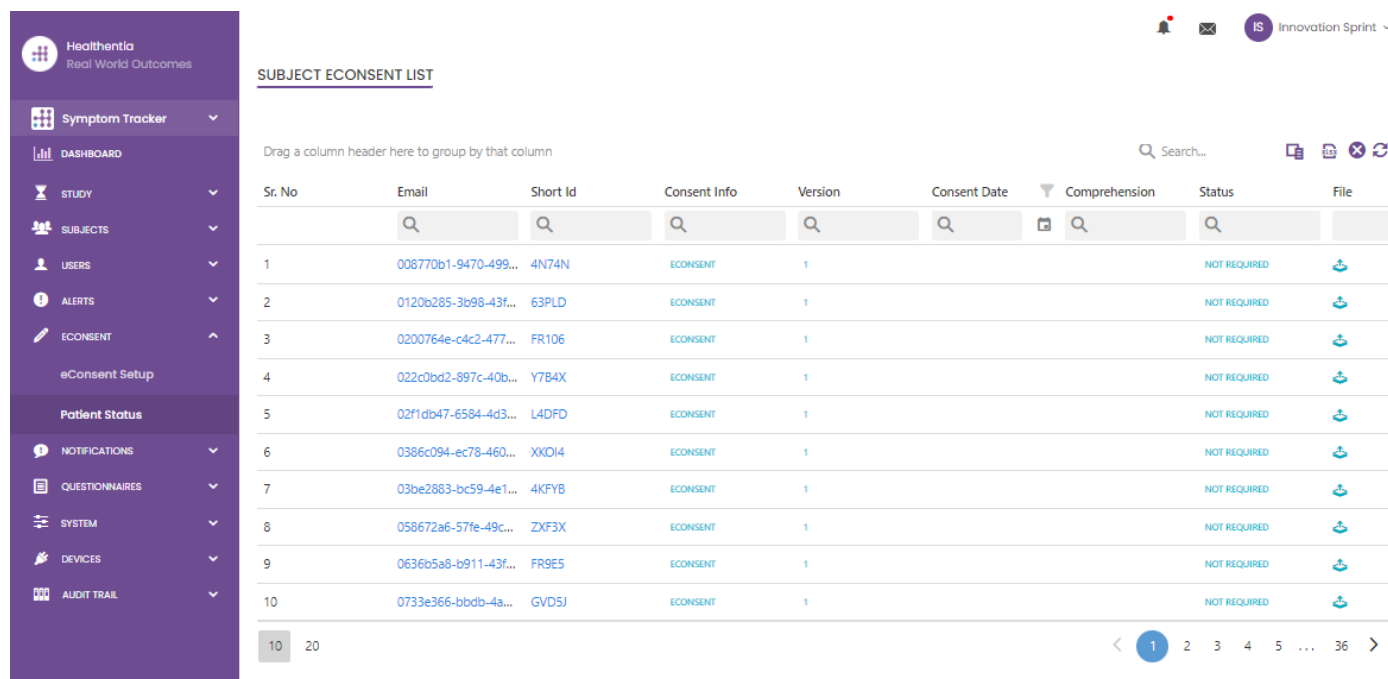
Title	Codename	Type	No. of Questions	Types	Created On	Status	Actions
Initial Questionnaire	ISPRINT_COVID19	Simple	11	Active	16/03/2020	ACTIVE	
Body temperature	ISPRINT_FEVER	Simple	1	ADD-DELETE-EDIT	18/03/2020	ACTIVE	
Cough	ISPRINT_COUGH	Simple	2	ADD-DELETE-EDIT	18/03/2020	ACTIVE	
Fatigue	ISPRINT_FATIGUE	Simple	1	ADD-DELETE-EDIT	18/03/2020	ACTIVE	
Difficulty of breath...	ISPRINT_DBREATH	Simple	1	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Headache	ISPRINT_HEADACHE	Simple	1	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Muscle Pain	ISPRINT_MUSCLES	Simple	2	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Abdominal pain	ISPRINT_ABDPAIN	Simple	1	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Diarrhea	ISPRINT_DIARRHEA	Simple	1	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Nausea	ISPRINT_NAUSEA	Simple	1	ADD-DELETE-EDIT	26/03/2020	ACTIVE	
Loss of Smell/Taste	ISPRINT_TASTESM	Simple	1	ADD-DELETE-EDIT	31/03/2020	ACTIVE	
Oxygen saturation	ISPRINT_OXYGEN	Simple	1	ADD-DELETE-EDIT	07/04/2020	ACTIVE	
Weekly Questionnaire...	ISPRINT_COVID19_D	Simple	4		08/04/2020	ACTIVE	
Add or take photo	PHOTO	Simple	1	ADD-DELETE-EDIT	24/04/2020	ACTIVE	
Emotional State	EMOSTATUS	Simple	1	ADD-DELETE-EDIT	06/05/2020	INACTIVE	

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The screenshot shows the 'ROLES & PERMISSIONS FOR ORG_1' page in the Healthentia web portal. On the left is a sidebar menu with options: Symptom Tracker, DASHBOARD, STUDY, SUBJECTS, USERS, ALERTS, ECONSENT, NOTIFICATIONS, QUESTIONNAIRES, SYSTEM, DEVICES, and AUDIT TRAIL. The main content area has a search bar and a table of roles. The table has columns for 'Name' and 'Actions'. The roles listed are: App Manager, Application, CRO Assistant, Investigator, Organization Admin, System Administrator, test, and Viewer. Each role has three action icons: a pencil (edit), a trash can (delete), and a plus sign (add).

Name	Actions
App Manager	[Edit] [Delete] [Add]
Application	[Edit] [Delete] [Add]
CRO Assistant	[Edit] [Delete] [Add]
Investigator	[Edit] [Delete] [Add]
Organization Admin	[Edit] [Delete] [Add]
System Administrator	[Edit] [Delete] [Add]
test	[Edit] [Delete] [Add]
Viewer	[Edit] [Delete] [Add]

The screenshot shows the 'CREATE ROLE' page in the Healthentia web portal. It features a 'Name' field with the placeholder 'The Role'. Below is a tree view for 'Role Access Permissions'. The tree starts with 'All' (selected) and branches into 'Dashboard', 'Export', 'Subject', 'Reported Events', 'Reported Questionnaires', 'Consent Actions', 'Instant Notification', 'Scheduled Notification', and 'Questionnaire'. Each branch contains a list of permissions with checkboxes. For example, under 'Subject', permissions include 'Create', 'Edit', 'Details', 'Personal Details', 'Delete', and 'View Email'. Under 'Reported Events', permissions include 'Export' and 'Details'. Under 'Reported Questionnaires', permissions include 'Export' and 'Notifications'. Under 'Consent Actions', permissions include 'Export'. Under 'Instant Notification', permissions include 'Create', 'Edit', 'Delete', 'Details', 'ViewContent', and 'Export'. Under 'Scheduled Notification', permissions include 'Create', 'Edit', 'Delete', 'Details', 'ViewContent', and 'Export'. Under 'Questionnaire', permissions include 'Create' and 'Edit'.

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The screenshot displays the Healthentia web portal interface. On the left is a purple sidebar with navigation options: Symptom Tracker, DASHBOARD, STUDY, SUBJECTS, USERS, User List, User Invitations, ALERTS, ECONSENT, NOTIFICATIONS, QUESTIONNAIRES, SYSTEM, DEVICES, and AUDIT TRAIL. The main content area is titled 'USER MANAGEMENT' and contains a table of users. The table has columns for Email Address, Role, Organization, Created On, and Actions. Below the table are pagination controls showing 5 and 10 items per page, and a page number 2. To the right of the table is a 'CREATE USER' form with fields for First Name, Last Name, Email Address, Password, Confirm Password, Role, Time Zone, and Sites. The form includes 'Cancel' and 'Submit' buttons at the bottom right.

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Search...

Email Address	Role	Organization	Created On	Actions
6afe7900-8d69-4030-b546-8c94f2165170@healthentia.com	Application	ORG_1	02/10/2020	[Edit] [More]
bb426815-dbe6-4182-862e-56139ef4d9cd@healthentia.com	System Administrator	ORG_1	01/03/2021	[Edit] [More]
ef9be18f-06ec-4bc1-b551-a8fada21ffe8@healthentia.com	System Administrator	ORG_1	02/10/2020	[Edit] [More]
f0c0b242-19b7-49b3-9b56-cb442ae32b95@healthentia.com	System Administrator	ORG_1	02/10/2020	[Edit] [More]
f8340709-2eb5-48d7-b755-38c5db756c91@healthentia.com	System Administrator	ORG_1	02/10/2020	[Edit] [More]

5 10 < 1 2 >

CREATE USER

First Name: * Last Name: *

Email Address: * Password: Confirm Password:

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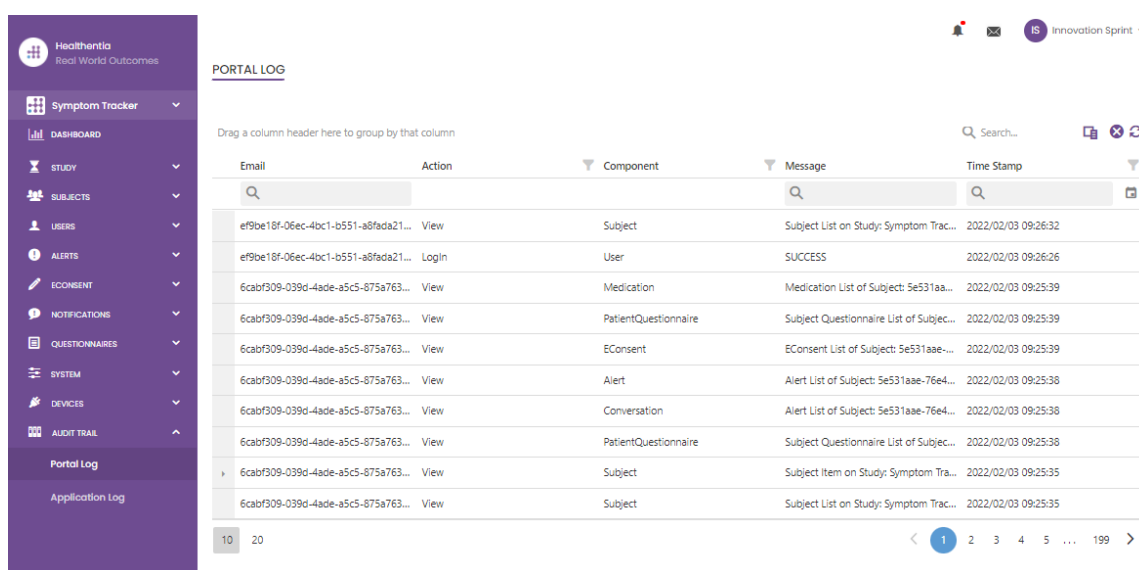
Email	Action	Component	Message	Time Stamp
ef9be18f-06ec-4bc1-b551-a8fada21...	View	Subject	Subject List on Study: Symptom Trac...	2022/02/03 09:26:32
ef9be18f-06ec-4bc1-b551-a8fada21...	Login	User	SUCCESS	2022/02/03 09:26:26
6cabf309-039d-4ade-a5c5-875a763...	View	Medication	Medication List of Subject: Se531aa...	2022/02/03 09:25:39
6cabf309-039d-4ade-a5c5-875a763...	View	PatientQuestionnaire	Subject Questionnaire List of Subjec...	2022/02/03 09:25:39
6cabf309-039d-4ade-a5c5-875a763...	View	EConsent	EConsent List of Subject: Se531aae...	2022/02/03 09:25:39
6cabf309-039d-4ade-a5c5-875a763...	View	Alert	Alert List of Subject: Se531aae-76e4...	2022/02/03 09:25:38
6cabf309-039d-4ade-a5c5-875a763...	View	Conversation	Alert List of Subject: Se531aae-76e4...	2022/02/03 09:25:38
6cabf309-039d-4ade-a5c5-875a763...	View	PatientQuestionnaire	Subject Questionnaire List of Subjec...	2022/02/03 09:25:38
6cabf309-039d-4ade-a5c5-875a763...	View	Subject	Subject Item on Study: Symptom Tra...	2022/02/03 09:25:35
6cabf309-039d-4ade-a5c5-875a763...	View	Subject	Subject List on Study: Symptom Trac...	2022/02/03 09:25:35

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